

The Acquisition SDK Programming Model allows Vantage *NXT* users to program Vantage *NXT* Research Ultrasound Systems without MATLAB[®] dependencies. The Acquisition SDK is comprised of a C-based API that offers familiar data structures and naming conventions, enabling easy porting of MATLAB scripts.


The Acquisition SDK is an ideal complement for users who aim to develop applications intended for commercialization or for integration of their legacy or third-party software to run on the Vantage *NXT* Platform.

The Acquisition SDK provides programmatic control of:

- Hardware communication
- Acquisition sequence definitions to acquire raw RF data only (no image reconstruction, processing, display or GUI included)
- Sequence loading, start and stop
- Direct Memory Access (DMA) management

And features:

- Hardware C-API with error checking for sequence correctness
- Example code
- Documentation

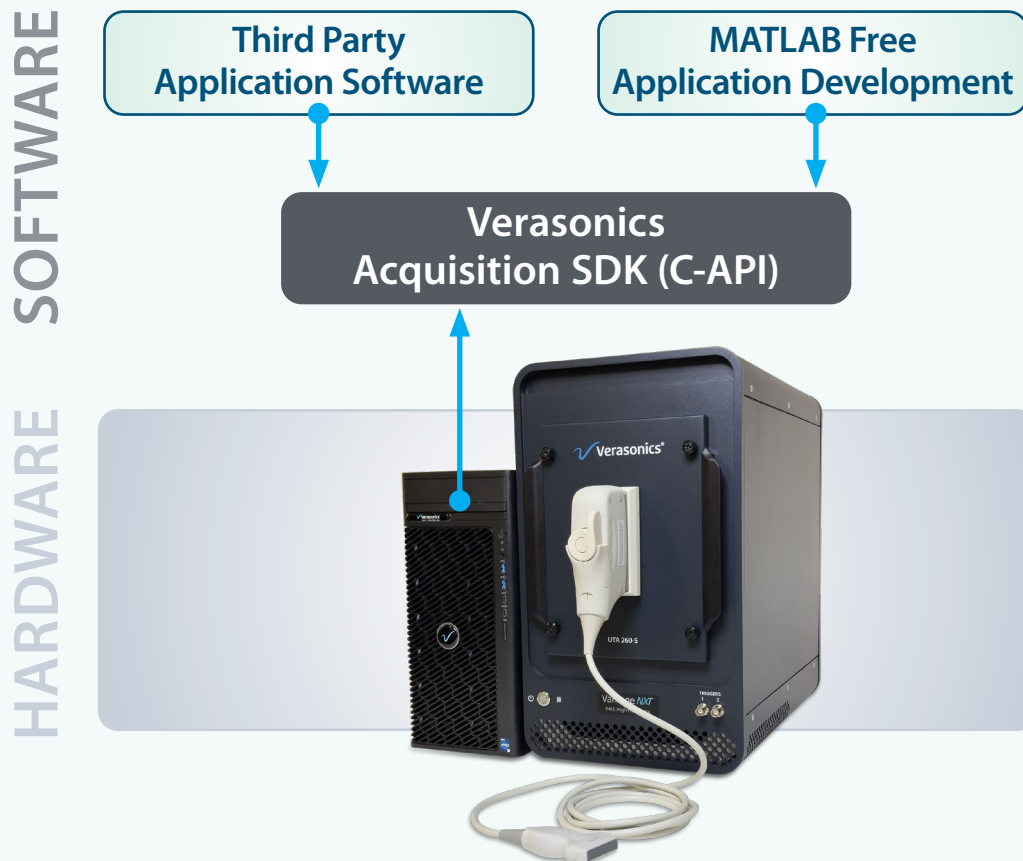

Supported functionality:

- Configurations supported: Vantage NXT 256, Vantage NXT 128
- Programming languages supported: C/C++ (users may utilize other languages, including Python)*
- Vantage NXT UTAs supported: 260-S, 260-D, 408-GE, 256 Direct, 64 LEMO, 128 LEMO, 160-DH/32 LEMO, 160-SH/8 LEMO, 160-SI/8 LEMO
- Licensable options supported: Triggers, Analog Inputs, Extended Transmit, GPU Toolkit with GPU Direct (for Linux OS only), HIFU Configurations

* Wrapper functions for other languages must be generated by the end user.

Complementary Programming Models to Increase User Productivity & Performance

MATLAB Programming Model Benefits:

- Rapid prototyping:
 - Leverage the flexibility of Verasonics' native image reconstruction
 - Develop custom processing functions written in other programming languages
- Utilize Verasonics' comprehensive set of example scripts, VSX GUI display, built-in Help System, and several other Verasonics-developed utilities and toolboxes
- Utilize MATLAB's built-in tools and capabilities

Acquisition SDK Programming Model Benefits:

- Implement sequences and algorithms prototyped in MATLAB
- Optimize code to run more efficiently without MATLAB dependencies
- Develop software using other programming languages
- Commercialize MATLAB-free products

Expanded Compatibility and Functionality

Capabilities to be added to the Acquisition SDK in 2025 include:

- Additional Vantage NXT configurations
- Additional Vantage NXT UTAs
- Additional Vantage NXT features

Not Included in Initial Release:

- Image Reconstruction
- GUI Display or Controls

Verasonics Inc.

11335 NE 122nd Way, Suite 100, Kirkland, WA 98034
www.verasonics.com | sales@verasonics.com



Verasonics reserves the right to change specifications without notice.
© Verasonics, Inc. 2025. All rights reserved. 4001-2025-61